

Bagmatic NOVO

NEW – PRODUCT PREVIEW



Features

- Best choice for field applications (battery driven, with no need for electricity)
- 8 hours of independent work, up to 30 donations autonomy
- 3-dimensional mixing of the blood donation
- OLED display with easy menu navigation
- Maximum and minimum flow detection ensures smooth treatment of the product
- Visual indication of error / warning / donation status
- Wireless connection for data transfer
- USB plug for data collection
- SD card data storage
- Powerful PC software for configuration of the device, data collection and analysis

Bagmatic NOVO

Bagmatic NOVO is our new blood collection monitor serving the special needs of mobile blood donations. As a light weight device with outstanding long battery life, it is developed to take a predefined blood donation and mix the product in a three-dimensional way. For quality assurance, Bagmatic NOVO monitors the entire blood donation and registers all significant events that may occur during the process.

Data collection

- With SD card
- USB connection
- Wireless connection



- Various connection possibilities


Bagmatic NOVO

Advantages

- Small size and extremely light weight
- High portability for mobile use
- Increased efficiency with all experience from Bagmatic SL
- Changeable battery system, no 220 V necessary
- Easy to clean and maintain
- Outstanding stability and long life time
- Bi-directional communication



- Data storage on SD card

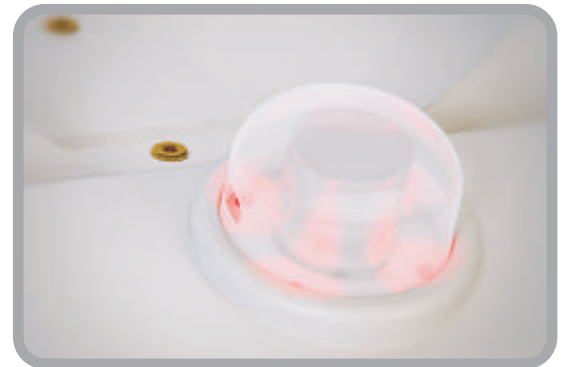

- Precise flow control



Technical data

Power supply:	14.4V Ni-MH battery pack
Power consumption:	standby: 80mA; during donation: 450mA
Battery capacity:	4000mAh
Battery charging time:	8h
Surrounding temperature:	10 – 50° C
Humidity:	maximum 85%
Dimensions:	300 x 220 x 140 mm
Weight:	3,0kg