


SEALMATIC

**NEW
product**

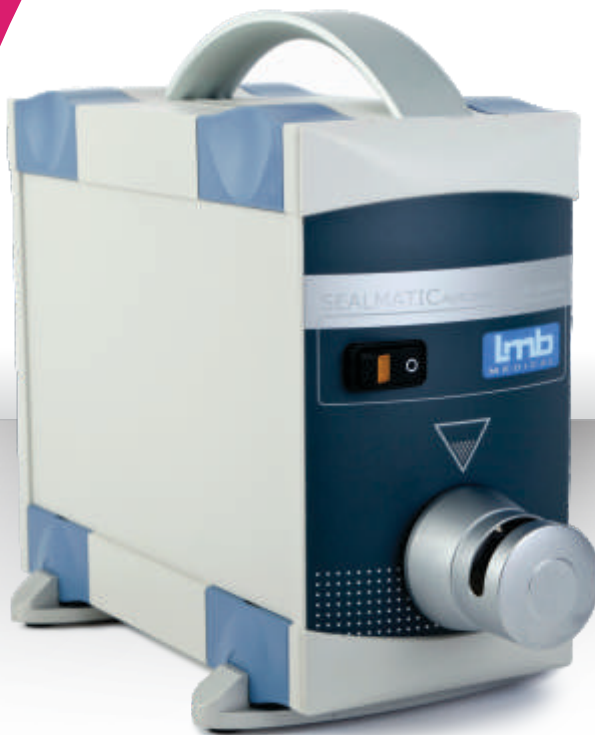


SEALMATIC is the latest answer to emerging demand of quality blood banking for high precision sealing combined with an outstanding short sealing time. The device SEALMATIC is designed as stationary sealer that offers just as well mobility in the donation room due to its low weight and comfortable carrying handle. Extremely low hemolysis figures complete SEALMATIC as a blood bank must-have.

- Mobile SEALMATIC for bloodbags
- Low hemolysis technology
- Shortest sealing time 0,99 sec



- High precision switch to activate the seal
- Easy to clean in case of blood leakage

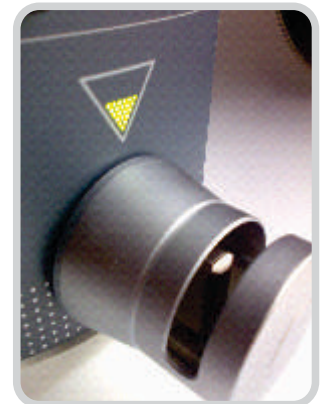


- Small head for transparency of the sealing process



- Precision electronics
- Precision mechanics
- Precision switching
- Precision hygiene



SEALMATIC

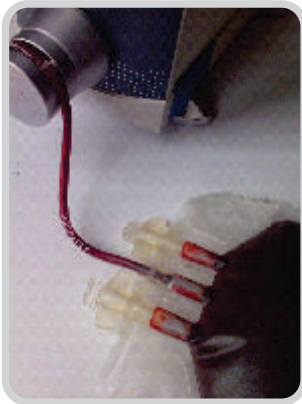


- Easy access for the user to insert the tubings
- Only 0.99 sec for one seal
- Tube diameter from 2,5 mm to 6,5 mm can be sealed

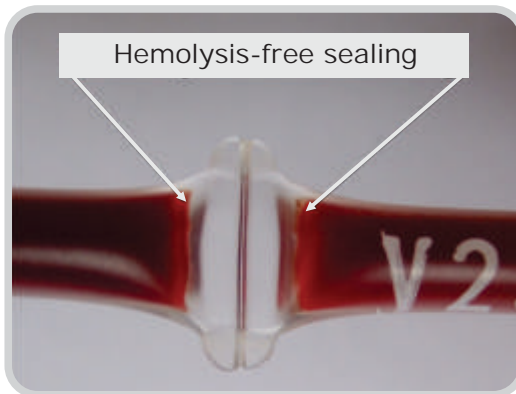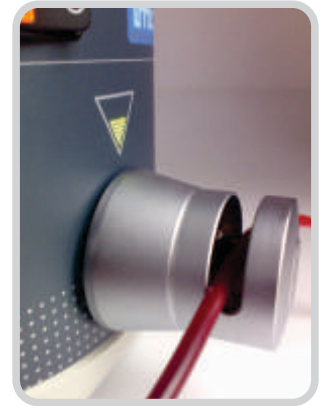

sample number	pre LDH-value before sealing in U/l	LDH-value after sealing in U/l
1	68	89
2	71	90
3	109	122
4	92	103
5	168	181
6	97	116
7	63	87
8	103	120
9	195	214
10	81	104
11	75	89
12	113	128
13	102	125



- High precision seal
- Optimal performance for separation

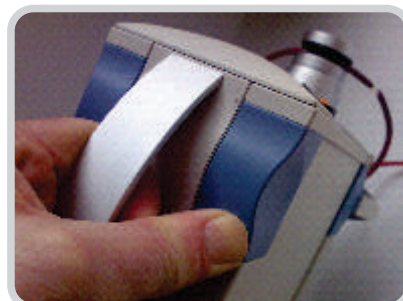

Specification

Sealing time	0,99 sec
Hemolysis	low hemolysis
Sealing power	50 W for 1 sec
Weight	3,3kg
Voltage	220V
Dimensions (mm)	w100 x h140 x d240